

PRESSCARE G300M

PRESSOTHERAPY | 2 OUTPUTS 3 INFLATION PROGRAMS

MEDICAL DEVICE SUITABLE FOR HEALTH, WELLNESS, BEAUTY AND RELAXATION NEEDS. INDICATED FOR LYMPHEDEMA, PROBLEMS OF VENOUS CIRCULATION, LEGS HEAVINESS, OEDEMAS, SWALLEN LIMBS, LYMPHATIC DRAINAGE, WATER RETENTION, CELLULITE, EMATOMS MESSAGES, VENOUS STASIS. IT CAN BE ALSO USED SIMULTANEOUSLY WITH THREE DIFFERENT ACCESSORIES. IT IS RECOMMENDED TO CONSULT THE DOCTOR BEFORE USING THE DEVICE.

TECHNICAL FEATURES

Outputs	2
Programs	3
Pressure	0-180 mmHg \pm 10%
Timer	up to 99 minutes
Mains supply	AC 230 V, 50 Hz
Size	185x180x230 mm
Weight	3 Kg

EQUIPMENT



- According to desired configuration
- 1 Power supply unit
- 1 Operating manual

PURCHASE CONFIGURATIONS

	BOOTS	SLEEVES	ABDOMINAL BAND	SPLITTING CABLE
· PRESSCARE G300M-1	1	/	/	/
· PRESSCARE G300M-2	2	/	/	/
· PRESSCARE G300M-3	2	/	1	1
· PRESSCARE G300M-1B	/	1	/	/



CE
0598

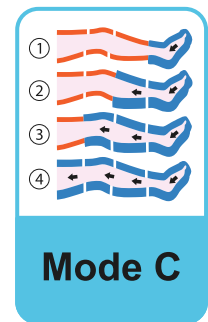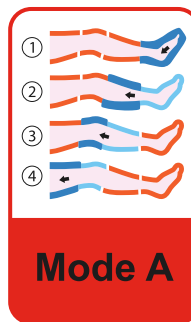

MODE A: according to the sequence of the program there will be an inflation and deflation of the boot one sector at a time, from the distal (foot or hand) toward the proximal part.

MODE B: according to the sequence of the program there will be an inflation of two chambers at a time, simulating the function of the overlap boots.

MODE C: the sequence of the program starts pressure from the feet and maintaining previous chambers inflated until the whole leg (or arm or abdomen) is inflated up to the groin.

PRESSCARE G300M



PRESSOTHERAPY | 2 OUTPUTS 3 INFLATION PROGRAMS

SIZE ACCESSORIES (4 SECTIONS)


Rev. En_01_2022

BOOTS SIZE "S"	(A) LENGTH GROIN - HEEL 72 cm (B) THIGH CIRCUMFERENCE 60 cm with EXTENSION 70 cm (C) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "M"	(A) LENGTH GROIN - HEEL 82 cm (B) THIGH CIRCUMFERENCE cm 67 con ESTENSIONE 77 cm (C) ANKLE CIRCUMFERENCE cm 37 con ESTENSIONE 42 cm
BOOTS SIZE "L"	(A) LENGTH GROIN - HEEL 92 cm (B) THIGH CIRCUMFERENCE cm 74 with EXTENSION 84 cm (C) ANKLE CIRCUMFERENCE cm 40 with EXTENSION 45 cm

SLEEVE	(A) 84 cm - (B) 63 cm (C) SHOULDER CIRCUMFERENCE 43 cm with EXTENSION 53 cm (D) WRIST CIRCUMFERENCE 31 cm with EXTENSION 36 cm
--------	--

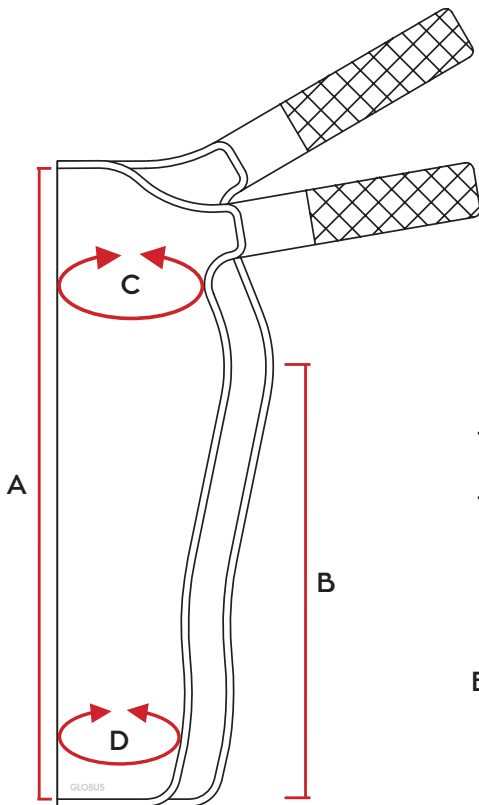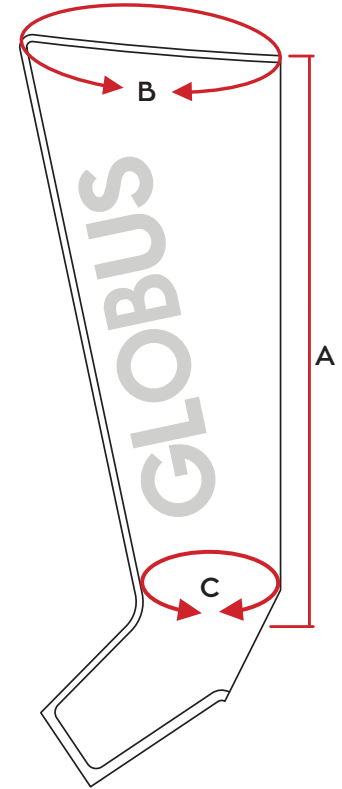

ABDOMINAL BAND	(A) 133 cm - (B) 38 cm ABDOMEN CIRCUMFERENCE from 82 cm to 114 cm
----------------	--

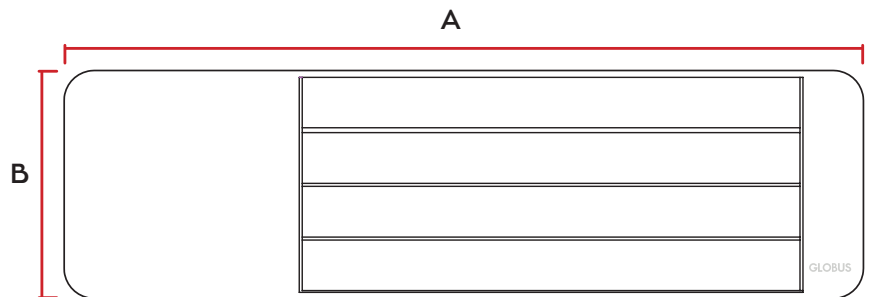

SPLITTER CABLE to use three accessories at the same time

SINGLE CABLE

N.B. All **EXTENSIONS** are to be considered an optional accessory available on request.

The **ACCESSORIES**, depending on the type of configuration, can be optional or included.